

Hip Stars

Journey's song "**Anyway You Want It**" is an apt description of **Sta-Rite™** Hip Guides.

<u>Listen</u>

Our Hip Stars offer Hardware and Pad options to meet any need:

- **Regular** provide maximum flexibility and function
- **Economy** provide strength and function at an affordable price
- **Bariatric** for those with extreme weight or tone needs.

Pads:

- 16 pads from Gel Ovations™
- Pads from other manufacturers

How to install? You won't believe it. You get a zipper bag that includes:

- ALL screws, T-nuts, and other hardware
- Instructions
- Everything for
 - power tilts from Quantum, Quickie, Motion Concepts, Merits, & Permobil
 - most popular manual tilts from Iris, Solara, Focus, & Fuze.

Remaining Problem is Pricing?

Don't let that get in your way. Call Chris **302-494-9070.** Or email Chris@GelOvations.com and we will figure out how to make this work for you and your clients.



Choosing the Right Hip Guide

Sta-Rite[™] Multi-Axis Adjustable Hip Guides come in 3 types so any client can be properly fitted.

The **most significant advantage** of the **Standard 6 Axis Adjustable Hip Guide** is it can be positioned <u>outside</u> <u>and over</u> the hardware so the patient is not squeezed or forced to buy a new, wider chair.

The **Standard 6 Axis Adjustable** & **Bariatric** have Track that allows them to be moved back and forth and to accommodate more complex client needs. The **Economy** version does not have a Track. To change the position of the **Economy** guide, it needs to be mounted in a different place along the length of the chair side and this is not always possible.

Aluminum cones are used for positioning in and out, and are **infinitely adjustable**. The interlocking steel crowns on the **Bariatric** adjust to fixed positions on the crowns. All the **Sta-Rite™ Hip Guides** can be easily adjusted **up** and **down**. Once positioned, the locking set clamp makes easy to return to the desired position.


The pads on **Sta-Rite™ Hip Guides** can be angled if desired. The Guides with Tracks can also be positioned forward or rearward along the leg. When forward, they can act as a Knee Button as well.

How to Pick Hip Guides or go to Gelovations.com and select How to Pick Hip Guides

Choosing the Right Sta-Rite™ Hip

| ТҮРЕ | Adjusts | Inward Adjustability including Pad | Outward Adjustability | Ability to Angle | Track | Material | Locking Type |
|---------------------------|---------|--|---------------------------|------------------------|-------|----------|--------------------------|
| Full Capability* SRMHG | 6 ways | up to 3.5" | out over the hardware* | Yes | Yes | aluminum | interlocking cones |
| Bariatric SRMHB | 6 ways | up to 2.5" | inward only | Yes | Yes | steel | interlocking crowns** |
| Economy SRMHE | 4 ways | up to 3" | inward only | Yes | No | aluminum | interlocking cones |

Sta-Rite™ Hip Guide codes: SRMHxyy

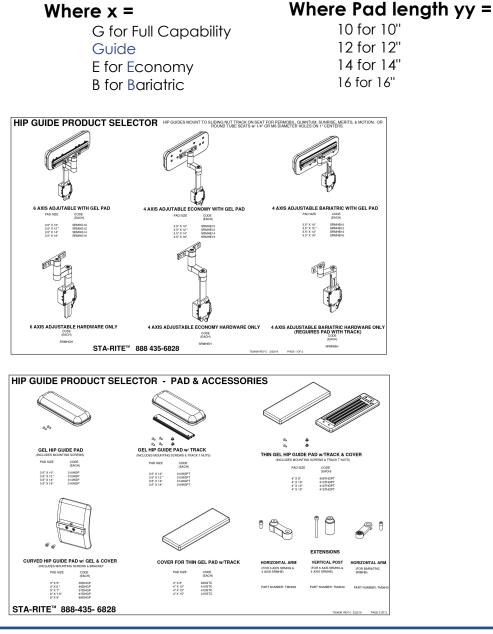



8245 Quebec Street Commerce City, CO 80022 Fax (888) 368-9598 Tel (303) 914-0141 Tel (888) 435-6828 service@GelOvations.com Falcon Rehabilitation

VA contract # 36F79718D0323

Expiration 12/31/27