

INSTRUCTIONS - SRASF1-03, SWING AWAY ADDUCTOR

1. ASSEMBLE & POSITION SWIVEL MOUNT TO SEAT TRACK.

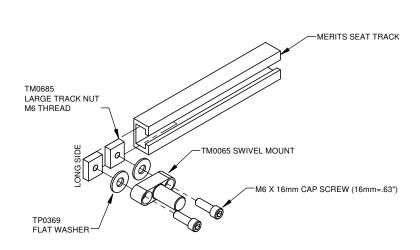

- 1. INSERT CAP SCREWS INTO SWIVEL MOUNT AND WASHERS AS SHOWN.
- 2. THREAD TRACK NUTS ONTO CAP SCREWS 1 TO 1 1/2 TURNS.
- 3. SLIDE TRACK NUTS INTO TRACK SLOTS AND POSITION SWIVEL CLAMP AS REQUIRED.
- 4. TIGHTEN CAP SCREWS TO SECURE SWIVEL MOUNT TO TRACK.

2. PLACE ADDUCTOR ASSEMBLY ONTO SWIVEL MOUNT AND ADJUST POSITION.

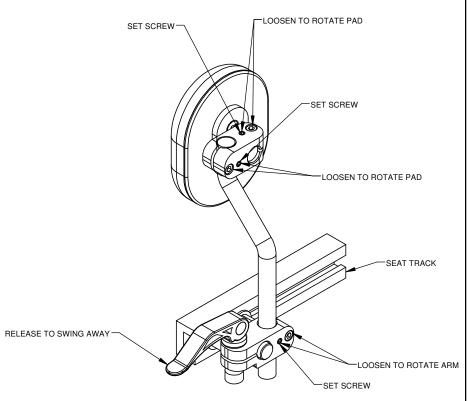

USE CAP SCREWS TO ADJUST POSITIONS.

TIGHTEN SET SCREWS (3 PLACES) AFTER FINAL POSITIONS ARE DETERMINED.

STA-RITE 888-435-6828

PAGE 1 OF 1